

DOCKAGE CHECKTEST

According to the Paperwork Reduction Act of 1995, no persons are required to respond to a collection of information unless it displays a valid OMB control number. The valid OMB control number for this information collection is 0580-0013. The time required to complete this information collection is estimated to average 20 minutes per response and 1 minute of recordkeeping, including the time for reviewing instructions, searching existing data sources, gathering and maintaining the data needed, and completing and reviewing the collection of information.

NOTE: TEST UNIT OPERATOR - FILL IN SHADED AREAS ONLY

DATE	FIELD OFFICE	AGENCY			LOCATION		
MACHINE SET UP	AIR NO.	FEED NO.	RIDDLE NO.	SIEVES	TOP NO.	MIDDLE NO.	BOTTOM NO.

TEST UNIT

MAKE		MODEL			SERIAL NO.			REMARKS
SAMPLE NO.	CLEARED		RIDDLE SEPARATION (Grams)	AIR SEPARATION (Grams)	THROUGH BOTTOM SIEVE (Grams)	TOTAL DOCKAGE (Grams)	OVER BOTTOM SIEVE (Grams)	
	HOPPER (Secs.)	SIEVE (Secs.)						
1								
2								
3								
AVG.								
TEST UNIT AVERAGE			%	%	%	%	%	
RESULTS BY:					DATE			


STANDARD UNIT

MAKE		MODEL			SERIAL NO.			REMARKS
SAMPLE NO.	CLEARED		RIDDLE SEPARATION (Grams)	AIR SEPARATION (Grams)	THROUGH BOTTOM SIEVE (Grams)	TOTAL DOCKAGE (Grams)	OVER BOTTOM SIEVE (Grams)	
	HOPPER (Secs.)	SIEVE (Secs.)						
1								
2								
3								
AVG.								
STANDARD UNIT AVERAGE			%	%	%	%	%	
RESULTS BY:					DATE			

DOCKAGE TESTER SUMMARY

	CLEARED		RIDDLE SEPARATION (%)	AIR SEPARATION (%)	THROUGH BOTTOM SIEVE (%)	TOTAL DOCKAGE (%)	OVER BOTTOM SIEVE (%)
	HOPPER (Secs.)	SIEVE (Secs.)					
TEST UNIT (Enter average results of test unit)							
STANDARD (Enter average results of standard unit)							
DEVIATION							
ALLOWABLE DEVIATION			0.10 %	0.10 %	0.10 %	0.15 %	
RESULTS							

REMARKS	
SUMMARY BY:	DATE

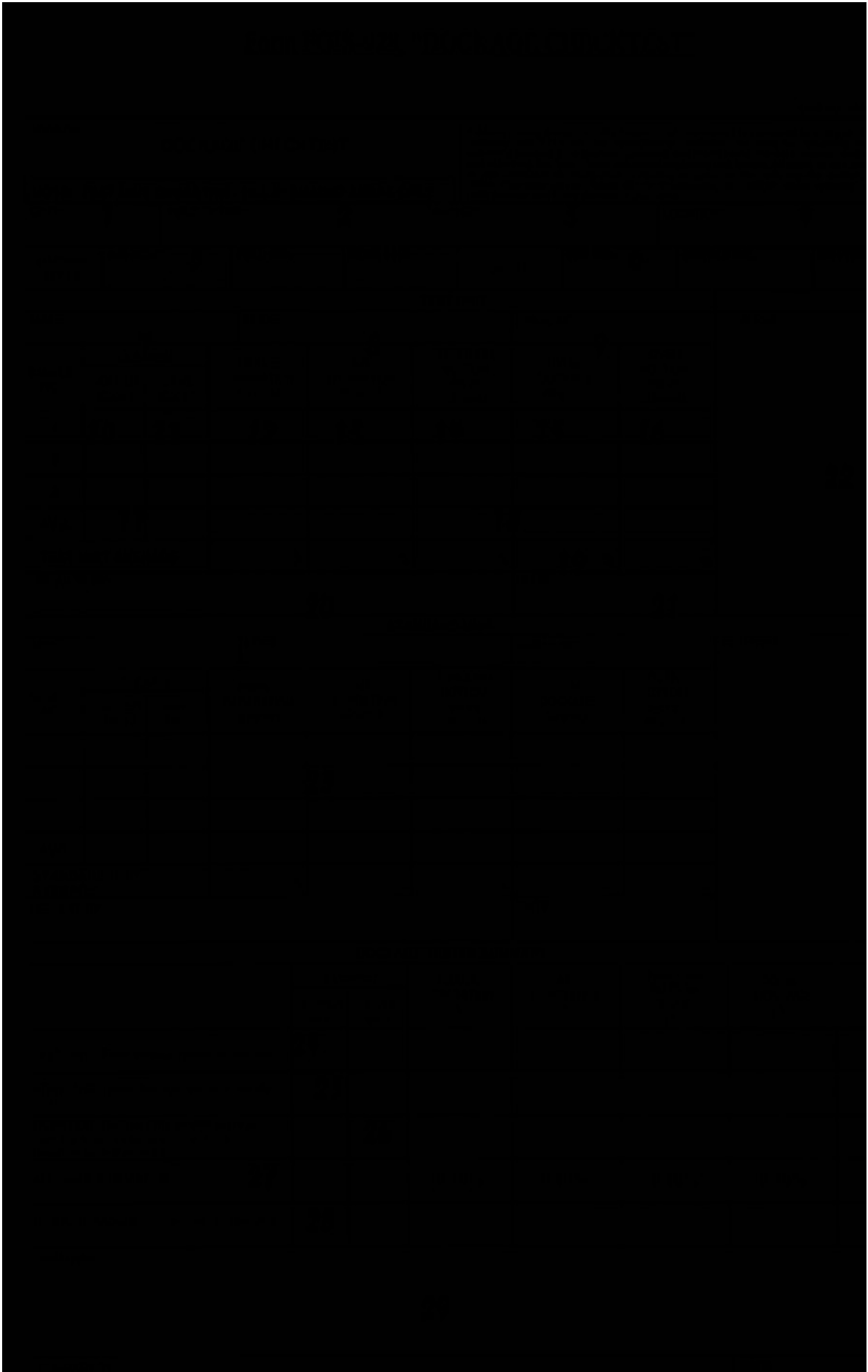



INSTRUCTIONS FOR COMPLETING
FORM FGIS-928, "DOCKAGE CHECKTEST"

- 1 Date mailed.
- 2 Field office that performed the test, when applicable.
- 3 Agency that performed the test, when applicable.
- 4 Location of test unit dockage tester.
- 5 For tests using wheat, show: air - 4, feed - 6, riddle - 2
For tests using corn, show: air - 1, feed - 10, riddle - n.
- 6 For tests using wheat, show: top - n/a, middle - 2 & I.I
For tests using corn, show: top - 3, middle - n/a, bottom
- 7 Make of the test unit dockage tester.
- 8 Model of the test unit dockage tester.
- 9 Serial number of the test unit dockage tester.
- 10 Time elapsed for each sample to clear the hopper, show
- 11 Time elapsed for each sample to clear the middle sieve
- 12 Weight of the separation that passes over the riddle, show
- 13 Weight of the separation that is removed by air, show
- 14 Weight of the separation that passes through the bottom
- 15 Total dockage (riddle, air, and through-the-bottom), show
- 16 Weight of the separation that passes over the bottom sieve
- 17 Averages of the test unit's hopper and sieve clearance
- 18 Averages of the test unit's test results, shown to 0.01 %
- 19 Convert the average separation weights to a percentage
- 20 Name of the person who tested the unit.
- 21 Date that testing of the test unit was completed.
- 22 Test unit operator's remarks and RPM information.
- 23 Standard unit's information and test results.
- 24 Test unit's average test results, shown to 0.01 %.
- 25 Standard unit's average test results, shown to 0.01 %.
- 26 Difference (24 minus 25). Show the plus or minus deviation
- 27 Allowable deviations.
- 28 Indicate IN or OUT of tolerance for each type of separation
- 29 Remarks
- 30 Name of person who completed the Summary section
- 31 Date the summary was completed.



.

/a.

D. No., bottom - 2 & I.D. No.

n - n/a.

wn in whole seconds.

z, shown in whole seconds.

own to 0.01 g.

1 to 0.01 g.

m sieve, shown to 0.01 g.

own to 0.01 g.

ieve, shown to 0.01 g.

times, shown in whole seconds.

z.

e of the total sample weight, shown to 0.001 %.

viation.

ation.

n.